


Digital Ultrasonic Cleaner (Sonicator)

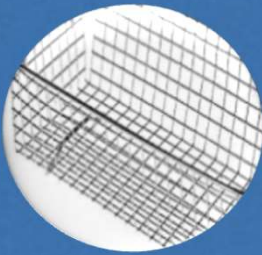

Product Description

The Digital Series is our most advanced ultrasonic cleaning unit. It features the benefits of heated cleaning (up to 69°C/156.2°F) plus complete led control with soft touch keypad capabilities and self-adaptive technology for consistent, worry-free operation. In addition to cleaning, baths can be used for a variety of applications. All sonicator cleaners are made with a stainless steel tank and chemically resistant plastic body. A cover is included with each unit.

A mesh basket or perforated tray must be ordered separately. A Basket is recommended for cleaning small items and the Tray is for use with heavier items or items that must be held upright and level (i.e. beakers).

Salient Features

- Inner Tank & Outer body made of Stainless Steel.
- Digital RED LED display.
- Temperature with timer function
- Ceramic heaters provides improved heating effect.
- Drainage and cool fan available in 6L to 25L Overheat protector, Cooling fan.
- Moisture-proof PCB.
- User friendly.
- High-performance transducers.

TECHNICAL SPECIFICATION

MODEL	BS-UDS3	BS-UDS6	BS-UDS9	BS-UDS12	BS-UDS25
Tank Capacity (Approx.)	2.5L	5L	8L	11L	24L
Tank Size (mm)	240x140x100	300x155x150	300x240x150	330x300x150	500x300x200
Overall Size (mm)	270x170x240	330x180x310	330x270x310	360x330x310	550x330x360
Ultrasonic Power	100W	150W	200W	300W	500W
Heating power	100W	300W	300W	400W	500W
Tank Drain	NA	Drain Valve			
Time Setting	1-99 min.				
Heating Temp.	0-80oC				
Display	Digital RED LED				
Frequency	40kHz				
Lid Material	SS 304				
Inner Tank Material	SS 304				
Outer Tank Material	SS 201				
Power Supply	AC 220V / 50 Hz				
Weight (approx.)	4 kg	8 kg	10 kg	13 kg	18 kg

STANDARD ACCESSORIES

- ✓Main instrument
- ✓Lid
- ✓Power Cord
- ✓Basket
- ✓Instrument

ORDERING DETAILS :-

SIZE, CAPACITY AND ORDERING DETAILS

MODEL	Tank Capacity (Approx.)	Tank Size (mm)	Overall Size (mm)	Ultrasonic Power	Heating power
BS-UDS3	2.5L	240x140x100	270x170x240	100W	100W
BS-UDS6	5L	300x155x150	330x180x310	150W	300W
BS-UDS9	8L	300x240x150	330x270x310	200W	300W
BS-UDS12	11L	330x300x150	360x330x310	300W	400W
BS-UDS25	24L	500x300x200	550x330x360	500W	500W



BIOGENIX SYSTEMS PRIVATE LIMITED

S-10, Pankaj Central Plaza, Pocket-5, Plot no-5,
Sector-12, Dwarka, New Delhi-110075, INDIA

Email: sales@biogenixsystems.com |

Phone: +91-11-43738899 | 9990251155-99 | 9911717057

Website: www.biogenixsystems.com

Dealer Details